


VascoPeel

VascoPeel is a radiopaque, pre-formable peel-away introducer set for use in the placement of pacing and defibrillator leads, balloon catheters, drainage catheters and other catheters through e.g. the vena subclavia.

VascoPeel is designed to be removed completely without displacement of the introduced catheter. The risk of kinking or traumatic effects at the insertion point is minimized by its optimized design.

The material ensures optimized pre-formability and stability which is important after removing the dilator from the sheath.

A white screw fixation connects the dilator with the sheath and ensures that the dilator remains fixed with the sheath when advancing the introducer into the vessel.

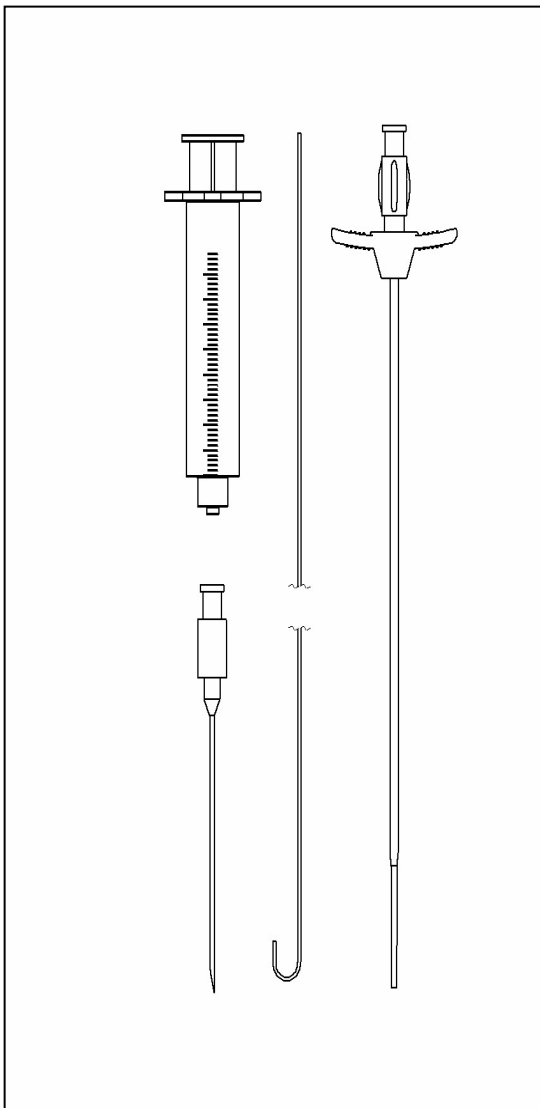

Technical Data:

Length (sheath)	14cm
Length (dilator)	17cm
French sizes	6F bis 12F

Package contents Set:

Set	Sheat/Dilator	Set with 2 Sheaths/Dilators
1 Peel-away Sheat introducer (14 cm)	1 Peel-away Sheat introducer (14 cm)	2 Peel-away Sheat introducer (14 cm)
1 Vessel dilator (17 cm)	1 Vessel dilator (17 cm)	2 Vessel dilator (17 cm)
1 0,035" x 45 cm safety J-tip guide wire		1 0,035" x 45 cm safety J-tip guidewire
1 thin-wall needle 18G		1 thin-wall needle 18G
1 LUER Slip syringe 10 ml		1 LUER Slip syringe 10 ml

Packaging Unit: 10 pcs. / Box

Artikel-Nummer:

Product	Set	Sheath/Dilator	Set with 2 Sheaths/Dilators
VascoPeel 6F	211028	211105	-
VascoPeel 7F	211091	211106	211549
VascoPeel 8F	211029	211283	211094
VascoPeel 9F	211030	211198	211095
VascoPeel 10F	211036	211199	211096
VascoPeel 10,5F	211155	211075	211097
VascoPeel 11F	211035	211200	211098
VascoPeel 12F	211156	211201	211099

The products are double packaged and sterile (ETO-Sterilization).